




# Product data sheet SintoGard

V 01 / 16.04.2014

# Wetrok SintoGard

## Disinfectants

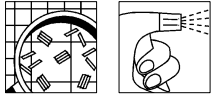

### Recommended dosage:



Contact time: see table,  
undiluted

### Benefits:

Ready-to-use, alcohol-based spray disinfectant.  
Evaporates without leaving any residues.  
Does not need to be rinsed off.

### Application:

For quick disinfection of small surfaces and objects (e.g. operating tables, examination chairs) and other alcohol-resistant materials.  
Can be used in the public healthcare sector, food industry and feedstuffs sector.

### Please note:

Treated surfaces must remain moist during the entire contact time.  
Only treat small surfaces (explosion hazards).  
Self-shining dispersions may be adversely affected by regular application. When using on paintwork or facing materials, it is advisable to test the product on a concealed surface beforehand.

Organism	Dosage	Time
Bactericide	Undiluted	30 seconds
Fungicide	Undiluted	30 seconds
Mycobactericide	Undiluted	60 seconds
Virucide (flu virus/influenza)	Undiluted	30 seconds

### Properties:

non-foaming

### pH value of concentrate:



### Important ingredients:

Ethanol 723 g/100 g

### Storage / shelf life:

Store in original container in a cool, well-ventilated place.  
Protect from frost. Keep container closed when not in use.  
Shelf life (unopened): 24 Months.

### Recommended application methods:

#### Surface disinfection

For well cleaned and dried surfaces (small surfaces and devices), spray SintoGard onto Quick & Clean Dry (cloth must be moistened thoroughly) and wet wipe entire surfaces/devices. Do not dry or rinse off. Wear mouth protection and gloves.